

visicolor® ECO Nickel

Níquel

470 nm


Method(e) / Método

5401


0.04–5.00 mg/L Ni²⁺

Test 5-40


REF 931 240

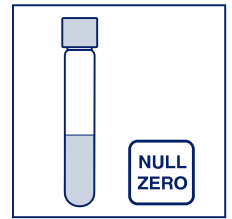


Probe / sample
échantillon / muestra

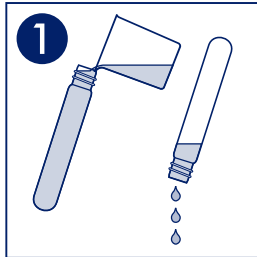


1
Küvetten spülen
Rinse test tubes
Rincer les éprouvettes
Lavar los tubos

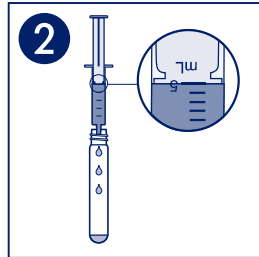


2
5 mL Probe
Sample
Echantillon
Muestra

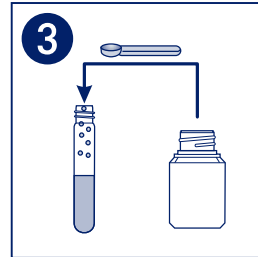

3
1 Ni-1

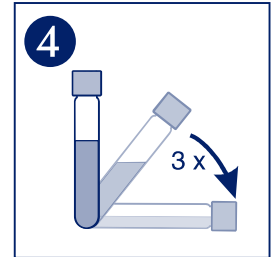


4
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

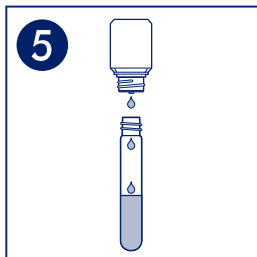


5
5 Ni-2



6
Schütteln
Shake
Agiter
Agitar

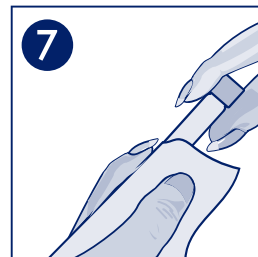


7
Säubern
Clean
Nettoyer
Limpiar

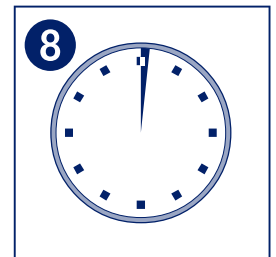

8
1'00 min



9
Messung
Measurement
Mesure
Medición

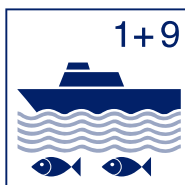

1+9
Meerwasser
Sea water
Eau de mer
Agua de mar