

MACHEREY-NAGEL
RNA, DNA, and
protein purification



NucleoSpin® RNA Set for NucleoZOL

- Benefit from both – NucleoZOL efficiency and NucleoSpin® purity
- Perfect recovery – no loss of total RNA, incl. miRNA
- Highest quality – even from difficult sample materials


NucleoSpin® RNA Set for NucleoZOL

The quick and safe solution for your NucleoZOL lysates

Product at-a-glance


Technology	Silica-membrane technology
Format	Mini spin columns
Sample material	Typically 500 µL NucleoZOL sample
Fragment size	Small RNA: 10–200 nt and large RNA: > 200 nt
Typical recovery	85–95 %
Typical RIN	Depending on sample quality
Elution volume	60 µL
Preparation time	< 1 h
Binding capacity	200 µg

Application data

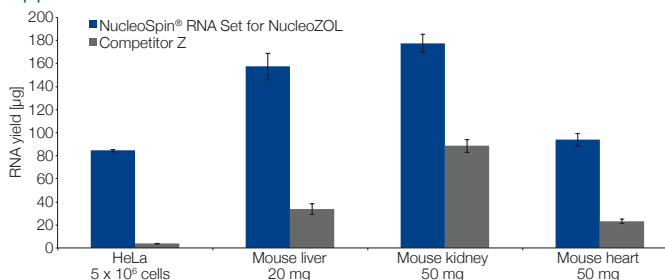

Total RNA was isolated with NucleoZOL combined with mini spin columns (NucleoSpin® RNA Set for NucleoZOL) or with a competitor product, that combines a liquid extraction with a column-based method (competitor Z).

Total RNA yields were quantified using an Agilent 2100 Bioanalyzer™.

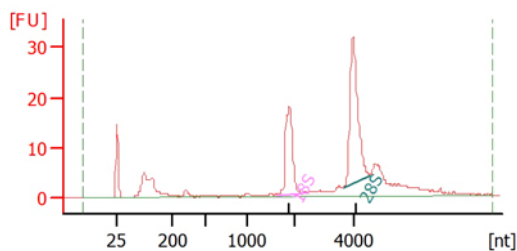

Total RNA was isolated with NucleoZOL in combination with the NucleoSpin® RNA Set for NucleoZOL from 5 x 10⁶ HeLa cells. RNA was analyzed on an Agilent Bioanalyzer™, a RIN of 9.8 indicates perfect RNA quality.

Ordering information

Product	Pack of	REF
NucleoSpin® RNA Set for NucleoZOL*	50 preps	740406.50
NucleoZOL	200 mL	740404.200

*Only be used in combination with NucleoZOL reagent; NucleoZOL not included.

www.mn-net.com

MACHERY-NAGEL



MACHERY-NAGEL GmbH & Co. KG · Neumann-Neander-Str. 6–8 · 52355 Düren · Germany

DE / International:

Tel.: +49 24 21 969-0
Fax: +49 24 21 969-199
E-mail: info@mn-net.com

CH:

Tel.: +41 62 388 55 00
Fax: +41 62 388 55 05
E-mail: sales-ch@mn-net.com

FR:

Tel.: +33 388 68 22 68
Fax: +33 388 51 76 88
E-mail: sales-fr@mn-net.com

US:

Tel.: +1 484 821 9984
Fax: +1 484 821 1272
E-mail: sales-us@mn-net.com

